


ENHANCE RDX

Enhance RDX Mattress System is a low to medium risk system for patients who are at risk of developing pressure ulcers according to your sound clinical judgment. Enhance RDX can also be used for patients who have an existing stage 1, 2, 3 and 4 pressure ulcers, in conjunction with your hospitals policy on pressure area management.

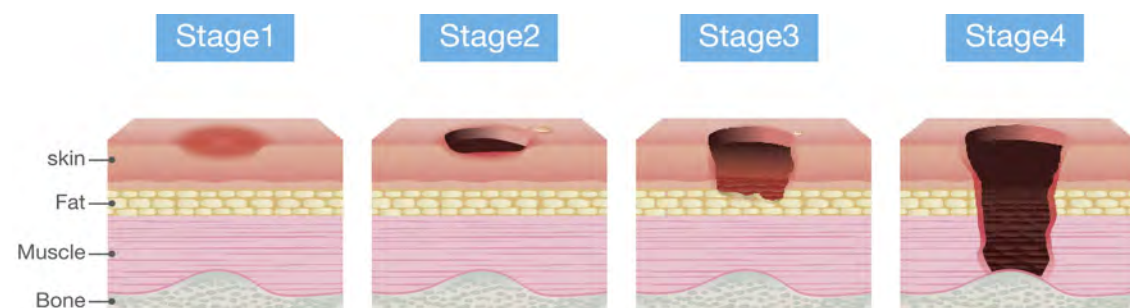

Optional Raised Bolster Cover available for added patient safety. Bolster cover has center access providing easy ingress and egress.

MASTER CONTROL UNIT	
Dimension:	(L)10.1"x (W)4.5"x (H)8.3"
Weight:	5.1 lbs
Rating Voltage:	110-120 volts
Max Current:	0.1A max.
Rating Frequency:	60Hz
Cycle Time:	10 mins
	Static function
	2 in 1 alternation
Min/Max Pressure:	12~47 mmHg +/- 5 mmHg
Flow-rate:	> 5 l/min
Protection Type:	Class I, Type BF
Low Pressure Alarm:	Visual / Audio
Mode Selection:	Static/Alternation Function

MATTRESS REPLACEMENT	
Model No:	Enhance RDX
Dimensions:	(L)80"x (W)36"x (H)8"
Weight:	20 lbs
Weight Capacity:	350 lbs
Cell Material:	Blue Nylon with PU coating
	Cell in cell design to stay inflated during power outage
Cell Number:	18 cells
Cover Material:	Blue Nylon with PU coating finish
Base Material:	Black Polyester with PVC coating
	CPR Quick Release System
HCPSC Code	E0277

PRIUS HEALTHCARE USA
Creating innovation and difference

T (813) 854-5464
E info@priushc.com
W priushcusa.com

160 SCARLET BLVD
OLDSMAR, FL 34677

ENHANCE RDX

Provides you, the caregiver, with peace of mind, and your patient with world class care and comfort




PRIUS HEALTHCARE USA
Creating innovation and difference

priushcusa.com

Perfect for patients and caregivers




Low slip base cover is a plus to keep the mattress in place on the bed

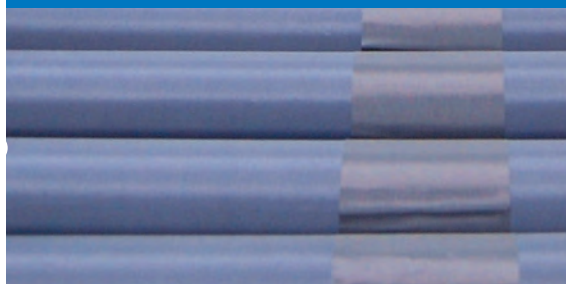

Welded vinyl bladder containment straps



Slimline pump design is unobstructive

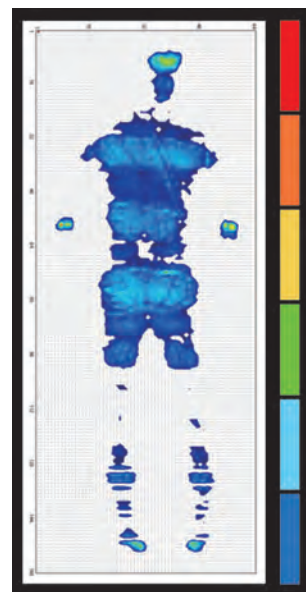

PRESSURE MAPPING

Enhance RDX provides constant low pressure by redistributing the patient's pressure over a greater area than a standard mattress.

ENHANCE RDX

- The Enhance RDX employs 8" air bladders with cell in cell design (4.5" of therapy, 3.5" of static air), to provide therapeutic alternating pressure and also allow the mattress to stay inflated in a power outage
- Durable woven polyester bottom cover with PVC lining for infection control
- The Enhance RDX replacement mattress can accommodate patient weights up to 350 lbs
- Low friction, low shear cover
- Lower foot end cells offload pressure in patient sensitive heel area



Enhance RDX is your partner in protection and treatment of pressure ulcers



Cable management keeps electrical cables off the floor



CPR quick release system deflates the system within 20 seconds

2-in-1 alternating air cycle ensures the maximum protection

The electronic control provides 8 adjustable levels of comfort depending on the patient's therapy requirement

The Power Failure Alarm is triggered during a power outage

The Low-Pressure Alarm is triggered when the air pressure of the mattress is lower than the default setup.



Hanging bracket fits on all type of beds

