



DUET


Duet Mattress System is a low to medium risk system for patients who are at risk of developing Pressure Ulcers according to your sound clinical judgment. Duet can also be used for patients who have an existing stage 1, 2, 3 and 4 pressure ulcers, in conjunction with your hospitals policy on pressure area management.

Stage1


Stage2

Stage3

Stage4

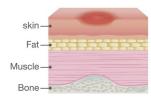



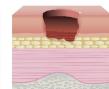







MASTER CONTROL UNIT	
Dimension:	(L)12.3"x (W)6.7"x (H)4.9"
Weight:	14.9 lbs
Rating Voltage:	110-120 volts
Max Current:	5A
Rating Frequency:	60Hz
Cycle Time:	3~95 mins
	Static / Pulsation function
	3 in 1 alternation / Lateral Rotation
Min/Max Pressure:	15~52 mmHg +/- 4 mmHg
	10 Digital Adjustable Comfort Control
Protection Type:	Class I, Type BF
Bottom Out Alarm:	Visual / Audio
Power Failure Alarm:	Visual / Audio
Features	User-Friendly LCD Control Panel
	Auto Fowler Detection
	Auto-Set Function
	Turning Function
	Patient Care Mode

MATTRESS REPLACEMENT	
Model No:	Duet
Dimensions:	(L)80"x (W)36", 42", 48"x (H)11"
Weight:	39.6 lbs
Weight Capacity:	700 lbs
Cell Material:	Clear TPU
Cell Number:	18 cells
Cover Material:	Blue Nylon with PU coating finish
Base Material:	Black Polyester with PVC coating
	CPR Quick Release System

T (813) 854-5464

E info@priushc.com

W priushcusa.com

160 SCARLET BLVD OLDSMAR, FL 34677



PRIUS HEALTHCARE USA


Perfect for



mattress base cover to reduce shearing and friction between mattress and bedframe



The modular construction of the cells makes for simplistic repair and cleaning

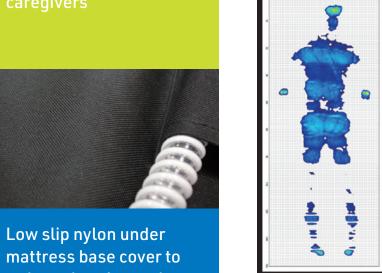


PRESSURE MAPPING

Duet provides constant low pressure by redistributing the patient's pressure over a greater area than a standard mattress.



- Equip with 2" convoluted foam base provides additional support during power failure situation
- Nylon PU cell design provides easy cleaning and minimizing the proliferation of bacteria
- Mattress replacement can accommodate patient weights up to 700 lbs

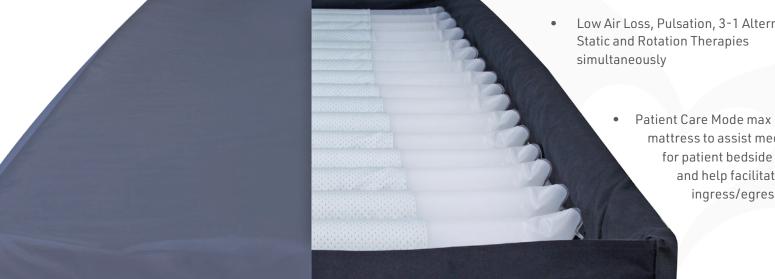

Low Air Loss, Pulsation, 3-1 Alternation,

 Patient Care Mode max Inflates mattress to assist medical staff for patient bedside treatment and help facilitate patient ingress/egress


Duet is your best choice for protection and treatment of pressure ulcers



Cable management keeps electrical cables off the floor



CPR quick release system deflates the system within 20 seconds

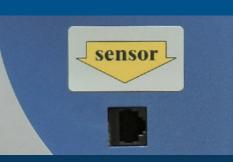

Incorporated with automatic sensing technology for fowler detection and pressure adjustment

Automatically sets optimal pressure based on patient's height and weight.

60 seconds of pulsation simulates a massage to assist in maximizing a patient's comfort.



Duet can detect changes in the angle of the bed frame during "Fowler Mode" and adjust the pressure based on these results.

Multiple lateral rotations provide alternative treatments to patients and medical staffs. Rotation therapy can be set from 3 to 95 minutes.



